


PRODUCT DATA SHEET

Gold NanoUrchins, >99% Endotoxin Free

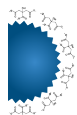

Citrate Capped

Catalog Numbers: UEF-50-XX, UEF-60-XX, UEF-70-XX, UEF-80-XX, UEF-90-XX, UEF-100-XX

Description

Cytodiagnosics Endotoxin Free Gold NanoUrchins have unique optical properties compared to spherical gold nanoparticles of the same core diameter. The spiky uneven surface causes a red shift in the surface plasmon peak and a larger enhancement of electromagnetic fields at the tips of the Gold NanoUrchin spikes compared to that of a spherical particle.

As an example, 100nm spherical gold nanoparticles have an SPR peak at 570nm while 100nm Gold NanoUrchins have a SPR peak at around 680nm.

In addition, binding of ligands such as proteins to the Gold NanoUrchin surface causes a larger shift in the surface plasmon peak compared to standard spherical gold nanoparticles.

The citrate-covered surface of our Gold NanoUrchins allows for efficient adsorption of primary antibodies and other proteins. In addition, Gold NanoUrchins can be further modified and functionalized through ligand-exchange with e.g., thiol-containing ligands such as PEG and oligonucleotides.

These particles can be used as an alternative to standard spherical gold nanoparticles in a wide range of applications such as electron microscopy, immunostaining and development of biological sensors.

Our Gold NanoUrchins are available in 6 different sizes and have uniform size distribution (CV <12%).

For custom sizes, formulations or bulk quantities please contact our customer service department.

Features



- Enhanced optical properties
- Endotoxin-free as tested by the Limulus ameocyte lysate (LAL) test for endotoxins and determined to be <0.5 EU/mL, whilst also being free from substances that interfere with the endotoxin testing reaction.
- Packaged with dust-free protection to maintain the highest level of sterility, also making them suitable for use in clean rooms
- Readily adsorbs proteins to the surface
- Ideal for cellular and cellular toxicity studies, immunological studies, experimental models, and aseptic assays

Applications

- Ideal for development of peptide and protein gold conjugates for use in applications such as blotting, lateral flow assays, LSPR assays, light microscopy, and transmission electron microscopy (TEM).
- Ideal for cellular and cellular toxicity studies

Characteristics

Core diameter: 50-100nm (Coefficient of Variance < 12%)
Polydispersity Index (PDI): < 0.20
Concentration: OD=1 (~ 0.05 mg/ml)
Absorbance (λ_{max}): 585-680nm (core diameter dependant)
Supplied in 0.1mM Phosphate-Buffered Saline (0.01X PBS)

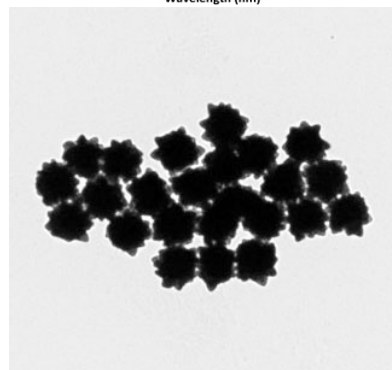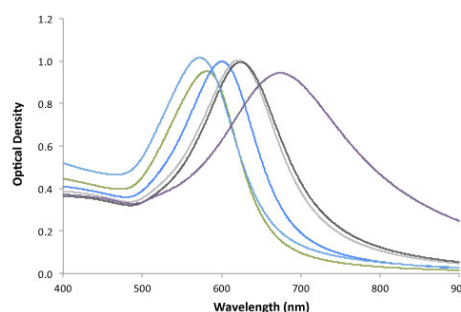

| Diameter (nm) | Peak SPR Wavelength (nm) | NPS/ml | Wt. Conc. (mg/ml) | Molar Ext ($M^{-1}cm^{-1}$) | Size Dispersity (+/-nm) | Particle Volume (nm ³) | Surface Area (nm ²) | Surface/Volume Ratio | Particle Mass (g) | Molar Mass (g/mol) | Molar Conc. |
|---------------|--------------------------|----------|-------------------|-------------------------------|-------------------------|------------------------------------|---------------------------------|----------------------|-------------------|--------------------|-------------|
| 50 | 585 | 3.51E+10 | 4.45E-02 | 1.72E+10 | <8% | 6.54E+04 | 7.85E+03 | 0.12 | 1.27E-15 | 7.64E+08 | 5.83E-11 |
| 60 | 585 | 1.96E+10 | 4.30E-02 | 3.07E+10 | <10% | 1.13E+05 | 1.13E+04 | 0.1 | 2.19E-15 | 1.32E+09 | 3.25E-11 |
| 70 | 600 | 1.20E+10 | 4.17E-02 | 5.03E+10 | <10% | 1.80E+05 | 1.54E+04 | 0.086 | 3.48E-15 | 2.10E+09 | 1.99E-11 |
| 80 | 620 | 7.82E+09 | 4.06E-02 | 7.70E+10 | <10% | 2.68E+05 | 2.01E+04 | 0.075 | 5.20E-15 | 3.13E+09 | 1.30E-11 |
| 90 | 630 | 5.37E+09 | 3.97E-02 | 1.12E+11 | <8% | 3.82E+05 | 2.54E+04 | 0.067 | 7.40E-15 | 4.46E+09 | 8.92E-12 |
| 100 | 680 | 3.84E+09 | 3.89E-02 | 1.57E+11 | <8% | 5.24E+05 | 3.14E+04 | 0.06 | 1.02E-14 | 6.11E+09 | 6.37E-12 |

* Data is approximated based upon a spherical nanoparticle. All concentrations at OD =1.

Protein Conjugation

A recommended starting protocol for conjugation of proteins to Cytodiagnosics standard gold nanoparticles can be found online at www.cytodiagnosics.com in the Technical Reference Section.

Storage

This product should be stored at 4°C. **DO NOT FREEZE.** If stored unopened and as specified, Cytodiagnosics Gold NanoUrchins are stable for at least 6 months.

Handling

When stored for a long period of time Gold NanoUrchins may sediment at the bottom of the flask, which is especially true for larger particle sizes. Prior to use, re-suspend the sedimented particles by swirling until a homogenous solution is obtained.

To maintain optimal performance, and stability of the colloidal gold, care should be taken to use clean storage containers if using other than supplied with the product.

Washing Gold Nanoparticles

Although it is not generally necessary to wash the Gold NanoUrchins prior to use, some applications may require additional washing procedures. The easiest way to remove possible contaminants in the nanoparticle solution is by centrifugation. Centrifugation force is dependent on size of the Gold NanoUrchins and should be adjusted according to Table I for optimal performance.

Note I: Since non-functionalized Gold NanoUrchins are sensitive to salt containing buffers, re-suspension should always be performed in ultra-pure water to prevent irreversible aggregation. Irreversible aggregation is characterized by a clear to bluish solution upon the addition of salt.

Note II: Please note that centrifugation can induce aggregation. To prevent aggregation, it may be necessary to add Tween 20 at a concentration of 0.025% w/v.

Procedure

1. Place 1ml aliquot of Gold NanoUrchins in a 1.5ml micro centrifuge tube.
2. Centrifuge the Gold NanoUrchins for 30 minutes using the appropriate G force determined by referencing Table I.
3. Remove the supernatant and re-suspend in an appropriate volume of ultra-pure water.
4. Vortex to re-disperse the particles.

| Size (nm) | Speed (g) | Time (min) |
|-----------|-----------|------------|
| 50 | 2,000 | 30 |
| 60 | 1,125 | 30 |
| 80 | 600 | 30 |
| 100 | 400 | 30 |

Table I. Appropriate G forces for centrifugation of Gold NanoUrchins. Note that recommended conditions are for a volume of 1ml and centrifugation using a microcentrifuge.

Product Safety and Handling

This product is for R&D use only, not for drug, household, or other uses. Please review the material safety datasheet (MSDS) available online for proper safety and handling procedures.

Ordering Information

For ordering call 866-344-3954 or visit us online.

| Catalog Number | Description | Sizes |
|-----------------------|---------------------------------------|----------------------|
| UEF-50-XX | 50nm Endotoxin Free Gold NanoUrchins | 20ml, 100ml or 500ml |
| UEF-60-XX | 60nm Endotoxin Free Gold NanoUrchins | 20ml, 100ml or 500ml |
| UEF-70-XX | 70nm Endotoxin Free Gold NanoUrchins | 20ml, 100ml or 500ml |
| UEF-80-XX | 80nm Endotoxin Free Gold NanoUrchins | 20ml, 100ml or 500ml |
| UEF-90-XX | 90nm Endotoxin Free Gold NanoUrchins | 20ml, 100ml or 500ml |
| UEF-100-XX | 100nm Endotoxin Free Gold NanoUrchins | 20ml, 100ml or 500ml |

*Indicates quantity, e.g., UEF-50-20 for 20ml of 50nm Endotoxin Free Gold NanoUrchins

For custom sizes, bulk quantities, and custom Gold NanoUrchin surface chemistry, please contact our customer service department.

Canada, Europe, Asia, Pacific and Africa

919 Fraser Drive, Unit 11, Burlington, ON Canada L7L 4X8
Tel: 866-344-3954 Fax: 289-204-9100
www.cytodiagnosics.com

United States, Mexico, South and Central America

5867 South Garnett Road, Tulsa, Oklahoma 74146 USA
Tel: 866-344-3954 Fax: 289-204-9100
www.cytodiagnosics-us.com